

Ambu[®] Man Compression



Key features

- Compact and portable training manikin
- Realistic training of chest compressions
- Visualisation of the correct compression depth
- Anatomical correct closed torso
- Easy placement of AED pads
- 5 pack for easy class training

AmbuMan Compression

Ambu Man Compression is a first aid and training manikin which offers human anatomical correct landmarks dedicated to the training of external chest compression.

The manikin allows an optimal training of the external chest compression, which is in the focus of the education of modern BLS algorithms.

Additionally the manikin can be used for realistic training of the placement of AED pads in order to get used in the handling of AED's.

The AmbuMan Compression 5 pack is optimal for class training. Optional the manikin can be equipped with head to give a more lifelike appearance.

Realistic chest compressions

With compression only CPR it is important to train the correct procedure for performing chest compressions. The manikin includes a feedback indicator that shows the correct compression depth.

Compact and portable

The AmbuMan Compression is a compact and portable training manikin.

The manikins are packed in a carrying of 5 manikins so that it can easily be transported.

Easy AED training

Due to the closed shaped torso, defibrillation pads can be placed and AED training can be performed realistically.

Specifications

Weight Torso	2,0 kg
Weight 5 x Torso including carrying bag	10,5 kg
Length Torso	45 cm
Max. compression depth	7,5 cm

Material


Base unit	ABS (Acrylnitril-Butadien-Styrol)
Chest skin	SEBS (Styrol-Ethylen-Butylen-Styrol)

Order information


Item no.	Description
A268 004 000	AmbuMan compression 5 pack



Closed torso

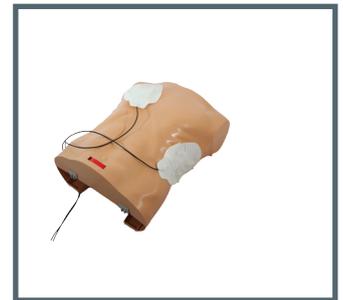

Placement of AED pads

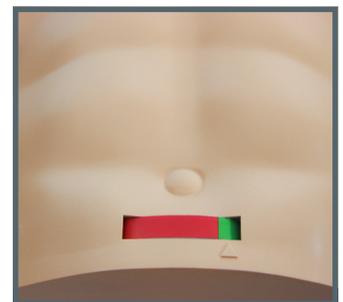

Compression depth indicator



Carrying bag for AmbuMan
Compression 5 pack

Ambu A/S
Baltorpbakken 13
DK - 2750 Ballerup
Denmark
Tel. +45 7225 2000
Fax +45 7225 2053
www.ambu.com